



13018 Cronos

Multi level stretcher with back rest and Trendelenburg




FEATURES

CRONOS is built in electropolished stainless steel, with foldable bed and side rails. The backrest and Trendelenburg are easily adjustable by means of gas springs, enabling virtually infinite positions. The 4 gas springs positioned below the loading bed allow perfect alignment at any height. Lastly, 4 wheels (Ø 200mm), 2 rotating wheels with brakes and 2 fixed swivel wheels, making CRONOS a handy and versatile stretcher.

MAIN FEATURES

Size	Adults
Usage	Transport
Multilevel	
Lying surface	With Trendelenburg
Wheel type	2 Rotating Wheels + 2 Big Fixed Wheels
Min height (cm)	48
Max height (cm)	86
Width (cm)	59
Length (cm)	198
Weight (kg)	37
Loading capacity (kg)	170
Material	Stainless steel ABS Polyethylene

IMAGES

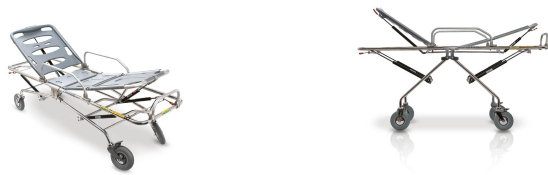

SEE ALSO



Icarus 13010
Multi level stretcher with side rails

ACCESSORIES



7013/N
3 pcs. Skay fire-proof electrowelded
mattress



7013/A
3 pcs. Skay fire-proof electrowelded
mattress



602/MEB
Bag with 3 belts 605/MEB
